CLINIC AIR

Instruction Manual Compressors SMART 2.25 SMART DRY 2.25 SMART 2.50 SMART DRY 2.50

C € EAE





INDEX

NTRODUCTION	4
STATEMENTS AND TESTING.	5
Statement of delivery of the instruction manual	5
Festing and commissioning certificate	6
EC DECLARATION OF CONFORMITY	
STATEMENT AND CE MARK	
Statement and CE mark	
GENERALI	
mportance of the manual	
Consultation notes	9 N
State "machine off"	n
Authorized operators	
Abbreviations	3
All rights reserved	3
Narranty14	4
TECHNICAL DESCRIPTION	
Name of models	
Names of parts	
Dimension's of the machine1	7
Technical data1	
ntended use of the machine1	
Source of energy supply1	
Products used for the operation	
Phases of operation2	0
Reasonably foreseeable misuse	
TRANSPORTATION AND HANDLING	
Iransportation of the machine	
Packaging2 Handling of packaging with a forklift or transpallet2	
Jnpacking	2
Vanual handling of the machine	2 2
Positioning of the machine	
INSTALLATION	
General warnings	
Pneumatic connection	
Electric connection	
SAFETY	
Safety devices	
Personal protective equipment (PPD)	
Residual risks	
Safety signs	Ó
JSE AND OPERATION	
Places occupied by authorized operators	
Pressure switch	2
Prior to ignition	3
Turning on the machine	3

J
•

Resettable thermal	34
Turning off the machine	35
Emergency stop	35
Activation after an emergency stop	36
Intervention of the safety valve	36
TROUBLESHOOTING	37
MAINTENANCE	38
Every week	39
Every month	40
Dryer filter ESC.F1	41
Dryer filter ESC.F3	42
Extraordinary maintenance	
DECOMMISSIONING	44
Storage	
Demolition	
SPARE PARTS	44

INTRODUCTION

Dear Customer, first of all, we wish to thank you for the confidence you gave us by purchasing your new "COMPRESSOR".

We are sure that your expectations will be met, given the level of technology achieved by our products with a continued commitment on a daily basis, which encourages us to grow, to know how to deal with the constant changes in technology expertise, manufacturing and sales.

We feel sure we can meet your future requirements, therefore, feel free to contact us. We will be happy to offer you all our experience and knowledge to answer any question you may have.

INFORMATION TO THE CUSTOMER			
DOCUMENT	PAGE TO DO		
STATEMENT OF DELIVERY OF THE INSTRUCTION MANUAL	HE 5 The authorised operators must com and sign to confirm the delivery and ding of the instruction manual.		
TESTING AND COMMISSIONING CERTIFICATE	6	Give to the test technician of the Manu- facturer or to the Authorized Dealer a copy of page 6 completed in all its parts.	
PRODUCT REGISTRATION FOR WARRANTY VALIDATION	14	See paragraph "WARRANTY" - Section 2) Page 14.	



STATEMENTS AND TESTING

STATEMENT OF DELIVERY OF THE INSTRUCTION MANUAL

The employer, to the extent required by applicable law on safety and health of workers at work, hereby declares it has delivered to the authorized operators this manual for the correct information and training on use and maintenance of the machine.

AUTHORIZED OPERATORS	DATE	SIGNATURE FOR RECEIPT
Mr		

Before filling in the table, photocopy it for future reference.



DATE				
Name:				
Street:				
City:				
State:				
This is to certify that the machine in question was tested by: Test technician of the manufacturer or authorized dealer:				
At the presence of the entrusted operator and Safety Representative: Mr				
Check for proper installation of the machine.				
g parameters.				
of all the safety devices.				
Commissioning and performance of a work cycle.				
Possible optimization of the working parameters.				
Machine stop.				
□ Delivery of the instruction manual.				
Safety Representative				



EC DECLARATION OF CONFORMITY

Gentilin S.r.l. - Via delle Tezze, 20/22 36070 Trissino (VI) - Italy VAT N.: 01262520248

Declares that the following machine:

SMART 2.25 - SMART DRY 2.25 SMART 2.50 - SMART DRY 2.50

Machine: Air volumetric compressor

Serial Number: _____

Complies with the following directives:

2006/42/EC Directive of the European Parliament and the Council of 17 May 2006 on machinery and amending Directive 95/16/EC (recast).

2014/30/UE Directive of the European Parliament and of the Council of 26 February 2014, replacing 2004/108/EC Directive of the European Parliament and of the Council of 15 December 2004, on the approximation of the laws of the Member States relating to electromagnetic compatibility and repealing Directive 89/336/EEC.

2014/35/UE Directive of the European Parliament and of the Council of 26 February 2014, replacing Directive 2006/95/EC of the European Parlament and Council, of December 12, 2006.

Person authorized to prepare the Technical Dossier: **Gentilin Giuseppe**, in quality of Technical Manager

Place: Via delle Tezze, 20/22 - 36070 Trissino (VI) ITALY

Date: _____

	Legal Representative (Gentilin Giampaolo) Signature:
C	entile las

Stamp:				
	1			 1.10
C	Re	ar	Ti	1)
	200			-



STATEMENT AND CE MARK


MARKING AND TECHNICAL DATA

The CE marking (Fig. 1) of the machine indicates compliance with the essential occupational health and safety rules set forth by the European Directives listed in the Declaration of CE Conformity.

The label is black and white printed on adhesive polymeric PVC and has the following dimensions: L = 70 mm - H = 55 mm.

It is applied externally on one side of the machine indicated in **Fig. 1** and states the following data, in a legible manner:

- LOGO, NAME AND ADDRESS OF MANUFACTURER
- CE MARKING
- DESIGNATION OF THE MACHINE (type / model)
- YEAR OF MANUFACTURE
- SERIAL NUMBER
- S1 VALUE
- MAX. WORKING PRESSURE (bar)
- VOLTAGE (V / Hz / A)
- SHAFT ROTATION SPEED (rpm)
- RATED POWER (kW)
- WEIGHT (kg)

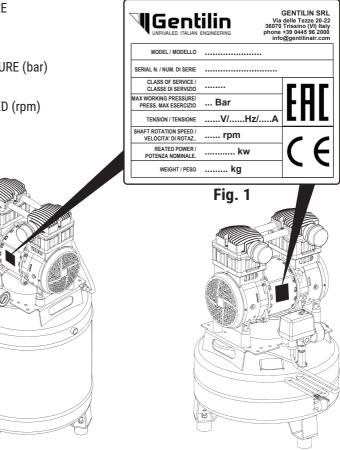

8

GENERAL

IMPORTANCE OF THE MANUAL

Before using the machine in question, the authorized operators are obliged to read and thoroughly understand this manual in all its parts.

This instruction manual has been prepared in accordance with the directions provided by the European Directives in order to ensure an easy and correct understanding of the topics covered by the authorized operators operating and maintaining the machine in question.

If, despite the attention being paid by the manufacturer while drafting this document, such operators should not fully understand this text while reading, they are requested to ask the manufacturer for the correct explanations and information in order to avoid incorrect personal interpretations that may jeopardize safety.

Before using the involved machine, the authorized operators are obliged to read and understand all parts of this manual and adhere strictly to the rules herein described in order to assure their own and others' safety, the best performance of the machine and ensure the maximum efficiency and durability of all its components.

This manual shall be available at any time for the authorized operators and shall be well protected and preserved, always close to the machine.

This manual should always be available for authorized operators and located in the vicinity of the machine well protected and preserved.

This manual must compulsorily be delivered together with the machine when it is transferred to another user.

The manufacturer disclaims all liability for injury to persons, animals and/or property caused by ignoring the rules and warnings described in this manual.

This manual reflects the state of the art of technique at the time of the machine marketing and cannot be considered improper just because, according to new experience, it may be updated later on.

In case of loss or damage of the manual request a copy to the manufacturer, specifying the identification data of the document: code, issuance and review.





CONSULTATION NOTES


SIMBOL	MEANING
	Signal of generic danger: It highlights risks to health and safety of operators and/or risk of damage or malfunction of the machine.
	Signal of generic obligation: It indicates a requirement (obligation to take action).
Θ	Generic prohibition sign: It highlights the prohibition to perform an action.
	Danger Sign EX: It highlights the risk due to explosion.

STATE "MACHINE OFF"

Before performing any type of maintenance and/or setting on the machine, disconnect the power supply ensuring that the switch (14) is in position "OFF (0)", the plug of the power cable (11) is unplugged and placed near the machine and the tank has been drained from residual compressed air by opening the drain valve (9).

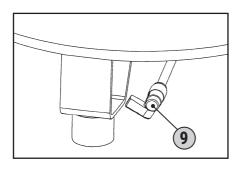

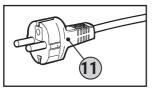

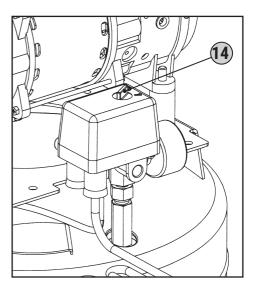

AUTHORIZED OPERATORS

D THE OPERATORS ARE AUTHORIZED TO MAKE ON THE MACHINE ONLY SPE-CIFIC ACTIONS OF THEIR COMPETENCE.

THE AUTHORIZED OPERATORS, BEFORE PERFORMING ANY WORK ON THE MA-CHINE MUST BE SURE TO BE IN POSSESSION OF FULL PSYCHO-PHYSICAL CA-PACITIES, SUCH AS TO ENSURE COMPLIANCE WITH THE SAFETY CONDITIONS.

This technical manual is intended exclusively for Authorized Operators to the use and maintenance of machines according to specific technical and professional skills required for the type of intervention. The symbols below are arranged at the beginning of a chapter and/or paragraph, to indicate which is the operator concerned with the subject matter.

OPERATOR IN CHARGE

It is a professionally trained operator, authorized to perform only the ignition, use, setup, tuning (with the mandatory protection enabled and the machine turned off) and the turning off of the machine and in total respect of the instructions in this manual, equipped with **personal protective equipment (PPD)** specified and occupying the places (see related paragraph).

OPERATOR IN CHARGE FOR HANDLING

It is a professionally trained operator, who has completed 18 years of age, in accordance with the law in the country of use, enabled to the conduct of forklifts, cranes or bridges wagons, to make safe transport, handling and unpacking the machine and/or parts of it, with the **personal protective equipment (PPD)** specified and occupying the places described (see related paragraph).

MECHANICAL/HYDRAULIC/PNEUMATIC MAINTENANCE TECHNICIAN

It is a qualified technician, authorized to perform actions only on the mechanical/ hydraulic/pneumatic parts to make adjustments, and maintenance and/or repairs, even with the protections disabled (upon consent of the Supervisor) in strict compliance with the instructions in this manual or other document specifically provided exclusively by the manufacturer, equipped with **personal protective equipment (PPD)** specified and occupying the places described (see related paragraph).



ELECTRICAL MAINTENANCE TECHNICIAN



It is a qualified technician (an electrician in possession of technical and professional qualifications required by the regulations) outboards to the cal equipment to make adjustments, and maintenance and/or repairs, even in the presence of voltage and with the protections disabled (upon consent of the Supervisor),

equipped with personal protection devices (PPD) envisioned and occupying the positions described (see related paragraph).

SAFETY REPRESENTATIVE



🔼 It is a qualified engineer appointed by the customer in possession of technical and professional gualifications required by the regulations on safety and health of workers in the workplace.

MANUFACTURER'S TECHNICIAN

🔊 It is a qualified technician made available by the manufacturer and/or authorized dealer to perform the required technical assistance, routine work and extraordinary maintenance and/or operations not mentioned in this manual that require specific

knowledge of the machine, equipped with personal protective equipment (PPD) specified (see related paragraph).



ABBREVIATIONS

Lists some abbreviations used in this book.

- cca. = About
- cap. = Chapter
- PPE = Personal protective equipment
- R = Right
- h = Hours
- EN = European Norm
- Es. = Example
- Fig. = Figure
- max. = Maximum
- min. = Minimum
- min = Minutes
- N. = Number
- page = Page
- par. = Paragraph
- Pos. = Position
- Rif. = Reference
- sec. = Seconds
- L = Left
- TAB. = Table
- see = SEE

ALL RIGHTS RESERVED

The reserved rights related to this instruction manual are in possession of the manufacturer. All information (text, drawings, diagrams, etc.) given here is confidential.

No part of this manual may be reproduced and distributed (in whole or in part) by any means of reproduction (photocopy, microfilm or otherwise) without the written consent of the Manufacturer.

All trademarks are property of their respective owners.



GUARANTEE

- 1) GUARANTEE DECLARATION: the Manufacturer gives an undertaking to the End-User to replace, repair or otherwise intervene in the machine, if it has conformity defects that compromise its correct use and operation, but only if such defects are attributable to the actual liability of the Manufacturer. The Manufacturer reserves the right to adopt the optimal solution to restore the conformity of the machine within a reasonable period of time. The guarantee does not cover parts of the Products subject to normal wear and tear and/or deterioration or defects resulting from inappropriate transport.
- 2) END-USER: the End-User is the person who buys the machine to use it directly. The End-User who buys the machine as a natural person for purposes that do not fall within the scope of commercial or professional activity is defined as a "consumer". The End-User is defined as a "professional or company" when the purchase is made by a professional or by a company acting for purposes which fall within the scope of its commercial or professional activity.
- **3) TERMINALS:** the Manufacturer is liable to the End-User when the lack of conformity occurs within the following terms from the date of purchase:

24 months if the purchase was made by a 'consumer' as described in point 2).

12 months if the purchase was made by a "professional or company" as described in point 2). All End-Users have the right to extend the guarantee (Extended Guarantee) for up to three years from the date of purchase of the machine if they: (a) have registered their purchase on the website www.gentilincompressors.com within seven days of purchase and (b) have completed the registration form completely and accurately, filling in the required fields with exact, complete and true information and have kept the registration code to be mentioned in the event of problems. Replaced parts or products will remain under guarantee until the original product guarantee expires.

- 4) VALIDITY: the guarantee is valid if the End-User notifies the Manufacturer of the conformity defect within 2 months from the date on which he has discovered it, but within the final term indicated in point 3) above. The above report must be accompanied by a regular purchase document (tax receipt or invoice). In the case of an Extended Guarantee, the declaration must also indicate the registration code.
- 5) FORFEIT: the guarantee will be void if the purchaser uses the machine improperly and/or not in accordance with the provisions of this "Instruction Manual" provided by the Manufacturer, or if the lack of conformity was caused by inexperience on the part of the purchaser or by chance. The use of the products after the discovery of defects releases the Manufacturer from any guarantee commitment or liability.
- 6) LIABILITY: the Manufacturer is released from any liability resulting from any damage sustained by the purchaser, resulting from loss or diminution of production which derives from any lack of conformity.
- 7) COSTS: the End-User has the right to the restoration of the conformity of the machine under guarantee without any costs for labour or materials. The shipping methods must be agreed with the authorised dealer or with the Manufacturer. In the event that the non-conformity is only apparent or the complaint arrives after the deadline, the Manufacturer may charge the cost of both transport and any repair that is indicated by the authorised dealer or the Manufacturer.



TECHNICAL DESCRIPTION

NAME OF MODELS

The machine is offered in 4 models. **Tab. 1** shows the model names and specifications that distinguish them.

MODEL	DESCRIPTION	
SMART 2.25	• Oil-free air compressor with 25 litres tank without dryer.	
SMART DRY 2.25	• Oil free air compressor with 25 litres tank with silica salt dryer.	
SMART 2.50	• Oil-free air compressor with 50 litres tank without dryer.	
SMART DRY 2.50	• Oil free air compressor with 50 litres tank with silica salt dryer.	


Tab. 1

NAMES OF PARTS

In figures are represented the major components that make up the machine.

KEY


- 1) Support feet
- 2) Non-return valve
- 3) Two-cylinder pumping unit
- 4) Electric motor
- 5) Handle
- 6) Air intake filters
- 7) Dryer/filter
- 8) Air tank 25 l
- 8a) Air tank 50 l
- 9) Condensate drain valve
- 10) Tank pressure manometer
- 11) Electric power supply cable with plug
- 12) Safety valve (calibrated at 10.5 bar)
- 13) Pressure switch
- 14) "ON (I) OFF (0) switch"

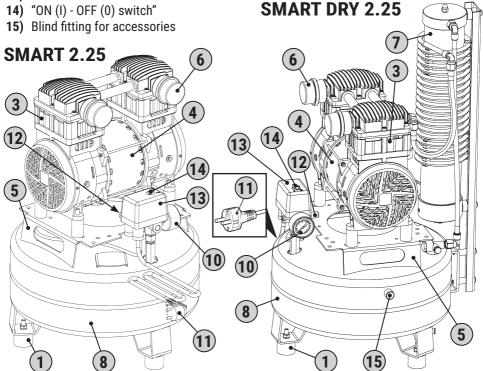



DIMENSIONS OF THE MACHINE

Tab.1 gives the clearance of the machine represented in Fig.1.

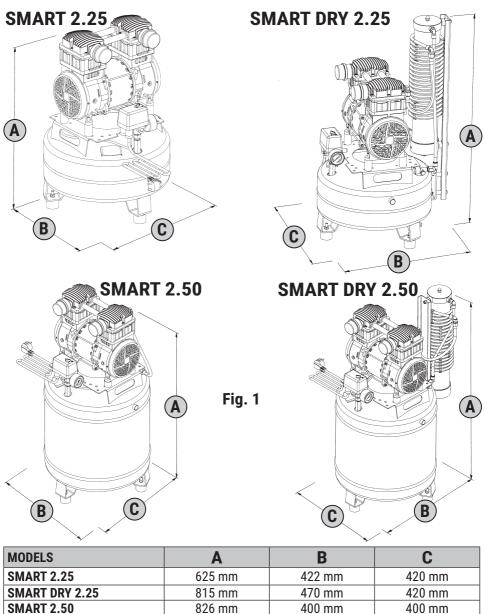

545 mm

826 mm

SMART DRY 2.50

400 mm

TECHNICAL DATA

MACHINE DATA

Tab. 2 shows the technical data which characterize the machine.

	SMART 2.25	SMART 2.50	
Power supply	230V/50 Hz		
Electrical Absorption	6,8 a	mper	
Units served		1	
Air delivered (5 bar) (liters/min)	140		
Intake air (5 bar) (liters/min)	235		
Tank capacity (liters)	24	48	
Nominal power (kW)	1,5		
Noise (dB(A))	67		
Weight (Kg)	32	??	
Ambient temperature (°C)	0-40		
Service class	70%		
Motor rotation direction	The compressor can be used in both directions of rotation		

SMART DRY 2.25	SMART DRY 2.50	
230V/50 Hz		
6,8 a	mper	
·	1	
14	40	
235		
24	48	
1,5		
67		
40	??	
0-40		
70%		
The compressor can be used in both directions of rotation		
	230V/ 6,8 a 14 23 24 24 1 6 40 0- 70 The compressor can be u	

Tab. 2



INTENDED USE OF THE MACHINE

The machine was designed and built for the use described in Tab. 1.

FIELD OF USE	Dentistry, laboratories or dental cabinets.
PLACE OF USE	Places indoors with sufficient lighting, ventilation, air environment having the characteristics referred to under par. "MACHINE DATA" and "PRODUCT USED FOR OPERATION" appropriate to the laws in force in the country of use for safety and health in the workplace. The machine must be placed on a plan that ensures stability in relation to the overall dimensions and weight, respecting the mi- nimum position shown in par. "POSITIONING OF THE MACHINE". WARNING! YOU MUST KEEP THE MACHINE OUT OF REACH OF UNAUTHORIZED PERSONS.
INTENDED USE	Compressed air production to operate dental tools.
OPERATORS IN CHARGE FOR USE	An authorized operator with the technical and professional capa- cities meeting the requirements.

Tab. 1

SOURCE OF ENERGY SUPPLY

Tab. 2 describes the source of energy supply of the machine.

ENERGY SOURCE	POWER SUPPLY	
ELECTRIC	Connection via plug on electrical switchboard. (See Technical Data)	

Tab. 2

PRODUCTS USED FOR THE OPERATION

Tab. 3 shows the product used for the operation.

PRODUCTS USED FOR THE OPERATION

Ambient air "clean" (good quality) free of pollutants harmful to health

Intake air filters (paper)

Dryer /filters

Tab. 3



PHASES OF OPERATION

In **Tab. 4** are listed the steps of operation of the machine (please see the relative paragraphs).

PHASE	OPERATION	
1	Turning on the machine.	
2	Adjusting the outlet pressure.	
3	If requested reset the resettable thermal.	
4	Turning off the machine.	

Tab. 4

REASONABLY FORESEEABLE MISUSE

The machine was designed and built for use provided, therefore it is absolutely prohibited any other type of application and use, in order to ensure at all times the safety of the Authorized operators and the efficiency of the machine itself.



RISK OF BURNS DUE TO ACCIDENTAL CONTACT WITH THE PUMPING GROUP AND ELECTRIC MOTOR AND PUMP UNIT. WARNING THERE IS A RESIDUAL RISK.



IT IS COMPULSORY TO INSTALL THE MACHINE IN A PLACE OF USE HAVING THE CHARACTERISTICS INDICATED ON PAR. "INTENDED USE OF THE MA-CHINE".

IT IS COMPULSORY TO USE THE MACHINE IN THE RESPECT OF THE LEGISLATION IN FORCE IN THE COUNTRY OF USE CONCERNING INTRODUCTIONS (NOISE).

IT IS COMPULSORY TO KEEP THE MACHINE OUT OF REACH OF UNAUTHORIZED PERSONS.

IT IS STRICTLY PROHIBITED THE USE OF THE MACHINE FOR MISUSE AND DIFFERENT USES FROM THE ONE INTENDED BY THE MANUFACTURER.

IT IS ABSOLUTELY FORBIDDEN TO USE THE MACHINE OUTDOORS.

IT IS STRICTLY PROHIBITED THE USE OF AIR TUBES (EXTENSIONS), FITTINGS AND NOT SUITABLE TOOLS AND/OR NOT CONFORM TO THE STANDARDS IN FORCE. IT IS ABSOLUTELY FORBIDDEN TO LIFT THE MACHINE WITH CRANE AND/OR FORKLIFT.

IT IS ABSOLUTELY FORBIDDEN TO LEAD AIR JETS AGAINST PERSONS AND/OR ANIMALS.

DO NOT CLIMB ONTO THE MACHINE.



IT IS STRICTLY PROHIBITED TO COMMISSION THE MACHINE IN ENVI-**RONMENTS WITH POTENTIALLY EXPLOSIVE ATMOSPHERE AND/OR IN THE** PRESENCE OF COMBUSTIBLE DUSTS (E.G.: WOOD POWDER, FLOUR, SU-GARS AND GRAINS).

IT IS ABSOLUTELY FORBIDDEN TO USE THE MACHINE IN AREAS CONTAINING SO-LID OR ABRASIVE PARTICLES.

TRANSPORTATION AND HANDLING

TRANSPORTATION OF THE MACHINE



💦 The machine is delivered to the customer through a "specialized transportation company" or by the manufacturer who, through its own personnel and appropriate means, in compliance with the regulations, shall ensure the operations of packaging, lifting, loading, transport and unloading in relation to the type of transport (land, sea or air).

ALL OPERATIONS OF LIFTING, LOADING, TRANSPORT AND UNLOADING OF THE MACHINE REQUIRE TO BE PERFORMED BY A "SPECIALIZED TRANS-PORT COMPANY" OR BY THE MANUFACTURER WITH PERSONNEL AND MEANS SUITABLE FOR USE.

PACKAGING



The machine is fixed to the pallet by straps and covered with a cardboard box. The packaging contains:



No. 1 compressor

No. 1 instruction manual

The characteristics of the packaging are indicated in Tab. 1.

Mod.	SMART 2.25	SMART DRY 2.25	SMART 2.50	SMART DRY 2.50
Packaging Weight	40 Kg	48 Kg	48 Kg	56 Kg
Pack measurement (machine + bed height)	68 × 57 × 95 cm			

Tab. 1

HANDI ING OF PACKAGING WITH A FORKI IFT OR TRANSPALL FT



🔊 IT IS MANDATORY THAT, IN THE RANGE OF HANDLING OPERATIONS. THERE ARE NO PEOPLE. ANIMALS AND/OR THINGS WHOSE SAFETY CAN BE COM-**PROMISED ACCIDENTALLY.**



THE HANDLING OF THE PACKAGE CAN BE PERFORMED ONLY BY THE USE OF A TRANSPALLET OR A FORKLIFT SUITABLE FOR USE. IT IS MANDATORY, DURING THE LIFTING OPERATIONS, (UP/DOWN), NOT TO MAKE SUDDEN MANOEUVRES AFFECTING STABILITY OF LOAD.


- 1) Insert the forks (2 Fig. 1) under the pallet (1), keeping them at the centre of the latter.
- 2) Raise the packaging (3) a few centimetres to make sure the load is actually centred in the forks.
- 3) Carefully move the package to the positioning area expected, therefore lowering the load until the full stability is reached.

UNPACKING

After laying the packaging on the ground do the unpacking as follows:

- 1) With a suitable tool, cut the straps (4) taking care not to be affected by the straps as a result.
- 2) Remove the cardboard box.

BE SURE TO DISPOSE THE PACKAGING MATERIAL ACCORDING TO THE DIFFE-RENT TYPES OF MATERIAL IN THE ABSOLUTE RESPECT OF THE LEGISLATION IN FORCE IN THE COUNTRY OF USE.

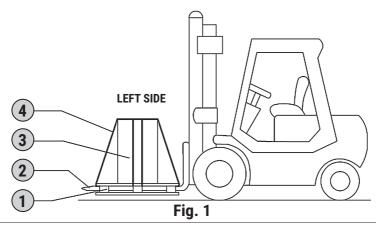



MANUAL HANDLING OF THE MACHINE

THE MANUAL HANDLING OF THE MACHINE MUST BE DONE BY AT LEAST TWO OPERATORS IN CHARGE, OR USING A SUITABLE MEAN (E.G.: TRAN-SPALLET), IN ACCORDANCE WITH THE REGULATIONS ON THE "MANUAL HANDLING OF LOADS", IN ORDER TO AVOID ADVERSE CONDITIONS THAT INVOL-VE RISK OF ERGONOMIC INJURIES TO THE BACK.

The machine can be moved manually by two operators in charge, using the handles (1 - Fig. 1) that are attached. Make sure the plug of the power supply cable is disconnected and that the cord is positioned above the tank.

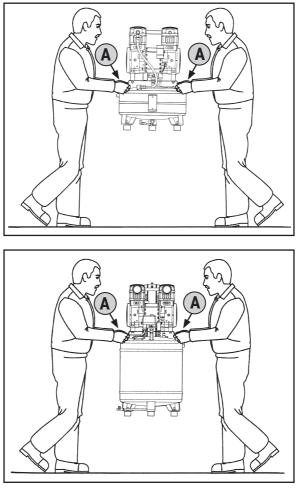


Fig. 1



POSITIONING OF THE MACHINE



The machine must be placed in a workplace with the characteristics described in par. "INTENDED USE OF THE MACHINE", placed on a flat floor that ensures stability in relation to the overall dimensions and to the weight.



IN ORDER TO ENSURE THAT THE AUTHORIZED OPERATORS OCCUPY THE WORK POSITIONS EXPECTED, IT IS MANDATORY TO ENSURE MINIMUM ME-ASURES OF POSITIONING AS SHOWN IN FIG. 2. IF THE ROOM TEMPERATU-RE, WHERE THE MACHINE HAS BEEN PLACED, EXCEEDS THE LIMITS SPECIFIED IN PAR."TECHNICAL SPECIFICATIONS", IT IS COMPULSORY TO INSTALL A SUITABLE **ADDITIONAL FAN, AS SHOWN IN FIG. 2.**

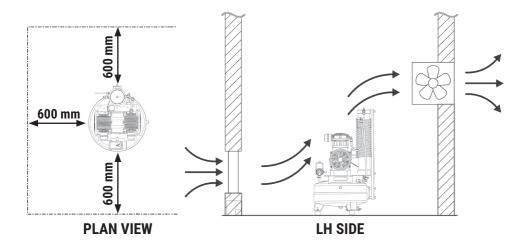


Fig. 2

INSTALLATION

GENERAL WARNINGS



IT IS COMPULSORY TO COMPLY WITH THE CORRECT PROCEDURES OF AS-SEMBLY/DISASSEMBLY AND CONNECTION OF COMPONENTS DESCRIBED IN THIS CHAPTER.

THE MANUAL HANDLING MUST BE IN COMPLIANCE WITH THE REGULATIONS ON THE "MANUAL HANDLING OF LOADS" IN ORDER TO AVOID ADVERSE CONDITIONS THAT INVOLVE RISK OF INJURIES TO THE BACK.

ALL OPERATIONS DESCRIBED BELOW ARE DONE TO PREPARE THE MACHINE READY FOR THE OPERATION.

PNEUMATIC CONNECTION



THE MANUFACTURER HAS NO LIABILITY FOR INJURY TO PERSONS, ANI-MALS AND/OR THINGS, RESULTING FROM FAILURE TO OBSERVE THE WAR-NINGS DESCRIBED ABOVE.



IT IS STRICTLY PROHIBITED THE USE OF AIR TUBES (EXTENSIONS), FIT-TINGS AND TOOLS NOT SUITABLE AND NOT CONFORM TO THE REGULA-TIONS IN FORCE.


IT IS COMPULSORY TO USE AIR TUBES (EXTENSIONS), FITTINGS AND TO-

OLS, ACCORDING TO WHAT IS FORECAST IN THE "INSTRUCTION MANUALS" PROVIDED BY THE MANUFACTURERS.



The machine must be connected upstream to a pneumatic system for the production sterilized air (by customer) connecting, to the air outlet **(Fig. 1)**, the tube for air inlet of the system equipped with a special fitting.



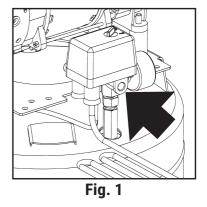

ELECTRIC CONNECTION

THE MANUFACTURER HAS NO LIABILITY FOR FAILURE OR MALFUNCTIONS OF OPERATION TO THE MACHINE CAUSED BY VOLTAGE DROPS EXCEEDING THE TOLERANCES PROVIDED BY THE DISTRIBUTING BODY (± 10% POWER - FREQUENCY ± 2%).

THE MAINS TO WHICH THE MACHINE IS CONNECTED MUST BE IN ACCOR-DANCE WITH THE REQUIREMENTS OF THE LAWS OF THE COUNTRY OF USE, MEET THE SPECIFICATIONS CONTAINED IN PAR. "MACHINE DATA" AND BE PROVIDED WITH A SUITABLE GROUNDING SYSTEM.

IT IS COMPULSORY TO INSTALL UPSTREAM OF THE POWER SUPPLY LINE AN ELECTRIC CUTTING DEVICE WITH A SUITABLE DISCONNECTING PROTECTION COORDINATED WITH THE GROUNDING PLANT.


ANY TYPE OF ELECTRIC MATERIAL USED FOR THE CONNECTION MUST BE SUITA-BLE FOR USE, MARKED "CE" IF SUBJECT TO LOW VOLTAGE DIRECTIVE 2006/95/ EC AND COMPLYING WITH THE REQUIREMENTS OF THE REGULATIONS IN FORCE IN THE COUNTRY OF USE OF THE MACHINE.

IF NECESSARY, IT IS COMPULSORY TO CONNECT THE MACHINE ONLY TO GENE-RATORS OF ELECTRICITY WITH A POWER GREATER THAN THE POWER SUPPLY INSTALLED, TO SUPPORT THE ABSORPTION CURRENT UPON STARTUP.

FAILURE TO FOLLOW THE WARNINGS ABOVE MAY CAUSE IRREVERSIBLE DAMAGE TO THE ELECTRIC EQUIPMENT OF THE MACHINE AND THE SUBSEQUENT VOID OF THE WARRANTY.



The machine can be connected to the mains by plugging the power cord **(11)** into the power outlet (by the Customer).

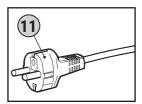

SAFETY

SAFETY DEVICES

The Safety Devices adopted are as follows (for the references on the drawings, please see the pages down):

REF.	DESCRIPTION	FUNCTION
1	Resettable thermal.	It is a thermal protection installed to the electrical sy- stem, which stops the electric motor in the instance of power overload and/or a short-circuit automatically
2	Tank pressure gauge.	It is a measuring device installed on the tank of the machine. It shows the air pressure inside the tank (0 \div 10 bar).
3	Safety valve.	It is a certified safety valve (set at 10.5 bar), installed in the pneumatic system under the pressure switch. It is used to discharge the overpressure of the system if the pressure switch, for any problems, does not function. The intervention of the safety valve forces the operator to turn off the machine and request the assistance of the maintenance personnel (see par. "INTERVENTION OF THE SAFETY VALVE").
4	Pressure switch.	It is an electro-pneumatic device (set at min. 6 bar, max. 8 bar) fitted onto the pneumatic plant. It is used to control the automatic start of the machine when the pressure drops to 6 bars and the automatic stop when the pressure reaches 8 bars.
5	Fixed cooling fan protective net.	Fixed protection made of plastic attached using screws fixed to the twin cylinder pump system base. It serves to stop accidental contact with the cooling fan when in motion.





IT IS COMPULSORY TO REPLACE ANY SAFETY DEVICES MALFUNCTIONING AND/ OR DAMAGED.


27

IT IS FORBIDDEN TO MANIPULATE, EXCLUDE, REMOVE AND/OR REPLACE ANY SAFETY DEVICE IN THIS MACHINE.

IT IS ABSOLUTELY FORBIDDEN TO REPLACE ANY SAFETY DEVICE OR COMPO-NENT WITH NON ORIGINAL SPARE PARTS.

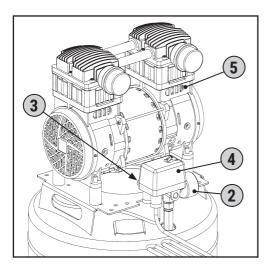

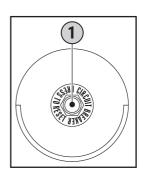

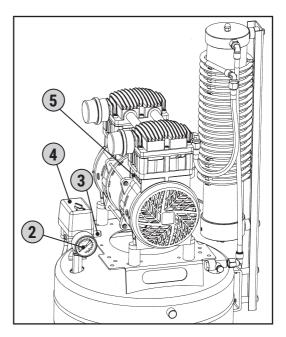



PERSONAL PROTECTIVE EQUIPMENT (PPD)

The use of personal protective equipment (PPD) is mandatory, in accordance with the legislation on safety and health in the workplace regulations of the country of use of the machine.

SEGN	PPE REQUIRED	TYPE OF APPLICATION
	Protection of hands (protective gloves against thermal hazards)	During handling, installation,
	Protection of feet (shoes with reinfor- ced toe)	operation and maintenance.


Tab. 1

IT IS COMPULSORY TO USE THE PPD SUGGESTED BY THE MANUFACTURER (SEE Tab. 1).

RESIDUAL RISKS

The authorized operator must remark that, although the manufacturer has taken every precaution technically possible to make the machine safe, in order not to compromise the functionality and productivity of the machine, there are still two potential residual risks described in **Tab. 2** and **3**.

RESIDUAL RISK N. 1	DANGER OF BURNS IN THE PUMPING GROUP AND ELECTRIC MOTOR	
SEVERITY OF THE INJURY	Light, normally reversible injuries.	
EXPOSURE	If the operator decides to perform a voluntary improper ac- tion, forbidden and not reasonably foreseeable.	
PROBABILITY	Low and accidental.	
STAGE WHERE THERE IS A RISK	During use and operation.	
MEASURES TAKEN	Safety signs. Compulsory use of personal protective equipment (PPD) and/or wait for the cooling of the machine.	

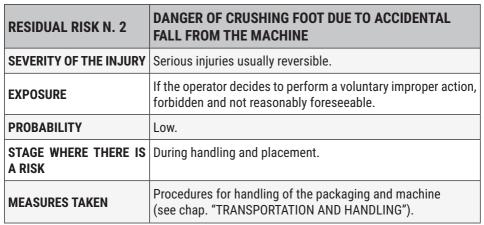

Tab. 3

SAFETY SIGNS

The Safety signs used are constituted by an adhesive label, applied externally to the machine (Fig. 1).



IT IS ABSOLUTELY FORBIDDEN TO REMOVE AND/OR DAMAGE THE SAFETY SIGNS APPLIED TO THE MACHINE.



IT IS COMPULSORY TO KEEP WELL CLEAN THE SIGNS IN ORDER TO ENSURE A GOOD VISIBILITY.

IT IS COMPULSORY TO REPLACE THE SAFETY SIGNS DETERIORATED, BY REQUE-STING THEM TO THE MANUFACTURER AND/OR AUTHORIZED DEALER.







SIGN	MEANING/PRESCRIPTION	SIGN	MEANING/PRESCRIPTION
4	DANGER: Electric Current		REQUIREMENT:
	DANGER: Automatic start		Read the instruction manual
			REQUIREMENT:
	DANGER: Dangerous Temperatures		Turn off the machine before carrying out maintenance


USE AND OPERATION

PLACES OCCUPIED BY AUTHORIZED OPERATORS

The Authorized operators, depending on the type of manual intervention required, must take the places shown in **Fig. 1**.



THE AUTHORIZED OPERATORS ACCORDING TO THE PLACE OCCUPIED, ARE OBLIGED TO PERFORM THE REQUIRED ACTIONS DESCRIBED IN THE RELA-TED SECTIONS.

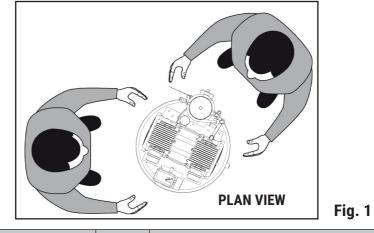

 AUTHORIZED OPERATORS
 PLACE
 DESCRIPTION OF INTERVENTION


 Image: Constraint of the supervisor of the supervisor, execute on the machine lifting operations, maintenance operations, handling and use of the machine.
 With the consent of the supervisor, execute on the machine lifting operations, maintenance operations, handling and use of the machine.




PRESSURE SWITCH

In **Fig. 1** is represented the control device to which the authorized operator must refer for the command and control of the machine.

REF.	NAME	FUNCTION	
1	Fitting.	Blind fitting for accessories.	
2	Tank pressure gauge.	It shows the outlet air pressure set through the regulator.	
3	Pressure switch. It is an electro-pneumatic device (set at min. 6 bar ma bar) fitted onto the pneumatic plant. It is used to contro automatic start of the machine when the pressure drops bar and the automatic stop when the pressure reaches 8		
4	Interruttore "ON (I)-OFF (0)"	Pos. "ON (I)": Turns on the machine; Pos. "OFF (0)": Turns off the machine.	
5	Safety valve.	It is a certified safety valve (set at 10.5 bar), installed in the pneumatic system under the pressure switch. It is used to discharge the overpressure in the system if the pressure switch, for any problems, does not function. The intervention of the safety valve forces the operator to turn off the machine and request the assistance of the main- tenance technicians.	

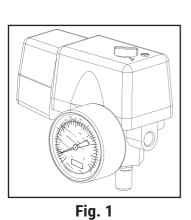

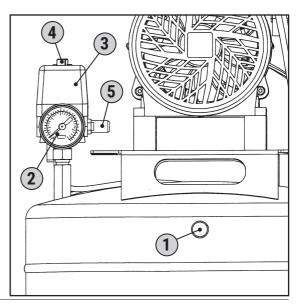



PRIOR TO IGNITION



BEFORE TURNING ON THE MACHINE THE OPERATOR IN CHARGE MUST MAKE THE CHECKS INDICATED IN TAB. 1.

No.	CONTROL	
1	Make sure that no unauthorized persons are near the machine.	
2	Make sure the machine is properly positioned.	
3	Make sure all connections are properly made.	
4	Ensure that the switch (14) is in pos. "OFF (0)".	
5	Make sure the condensate drain valve is closed (9).	
6	Use the Personal protective equipment (PPE) required.	
7	Make sure you have read and understood all parts of this "Instruction Manual".	
8	Make sure that after last use routine maintenance has been made.	

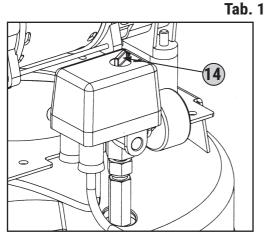

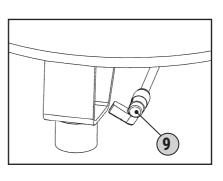

TURNING ON THE MACHINE



THE AUTHORIZED OPERATORS CAN TURN THE MACHINE ON ONLY AFTER THEY HAVE PERFORMED AS A COMPULSORY OPERATION. THE CONTROLS DESCRIBED IN THE CONTROLS TABLE BEFORE USING THE MACHINE IN **OUESTION IT IS COMPULSORY THAT THE AUTHORIZED OPERATORS READ** AND THOROUGHLY UNDERSTAND IN ALL PARTS THIS MANUAL. IT IS COM-PULSORY TO KEEP THE MACHINE OUT OF REACH OF UNAUTHORIZED PERSONS.



IT IS FORBIDDEN TO THE AUTHORIZED OPERATORS TO LEAVE THE MACHI-NE WITHOUT SUPERVISION DURING THE MAINTENANCE WORKS. THE MANUFACTURER HAS NO LIABILITY FOR INJURIES TO PERSONS, ANI-MALS AND/OR THINGS CAUSED BY THE INOBSERVANCE OF THE RULES AND IN-FORMATION DESCRIBED HEREIN.

Turn on the machine by moving the switch (14) in pos. "ON (I)" (the machine runs up to the max. operating pressure at 8 bar and then stops automatically).


RESETTABLE THERMAL


If there is an overload and/or a short circuit in the electrical system of the machine, the overload protection intervenes stopping the electric motor.

The reset will occur automatically when the temperature returns within the parameters preset by the manufacturer:

- 1) Set the switch (14) in pos. "OFF (0)".
- 2) Press the button of the resettable thermal (5).
- 3) Before restarting the machine wait a few minutes.

IF AFTER PERFORMING THE RESET THE MACHINE WILL NOT RESTART, THE OPERATOR IS REQUIRED TO OBTAIN THE INTERVENTION OF MAINTENANCE TECHNICIANS AND/OR THE AUTHORIZED DEALER.

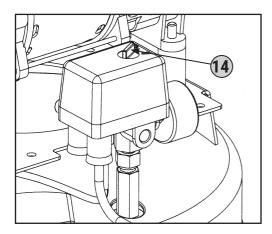

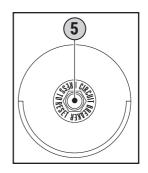

TURNING OFF THE MACHINE



At the end of the work cycle perform the turning off of the machine as follows:


- 1) Set the switch (14) in pos. "OFF (0)".
- 2) Unplug the power cord (11) from the socket and place it next to the machine.
- 3) Only in case of work cycle end perform the condensate discharge from the tank, acting on the condensate drain valve (9).


EMERGENCY STOP




The Emergency Stop of the machine can be done by bringing in pos. "OFF (0)" the switch (14). To avoid situations of imminent danger, the Authorized Operators are obliged to perform the following operations:

1	RUN IMMEDIATELY THE TURNING OFF OF THE MACHINE AS DESCRIBED IN PAR. "TURNING OFF THE MACHINE".
2	REPORT IMMEDIATELY THE EMERGENCY TO THE SUPERVISOR IN CHARGE.

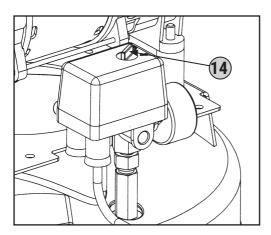

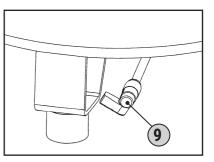

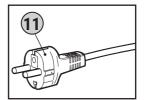





Only after removing the cause of the energency and here restore the ma-Only after removing the cause of the emergency and having carefully eva-🕥 the machine, with the consent of the Supervisor in charge restore the ma $igoplus_{1}$ chine as described in the following paragraph:

- 1) Set the switch (14) in pos. "OFF (0)".
- 2) If necessary, repeat the operations described in par. "TURNING ON THE MACHINE".

INTERVENTION ON THE SAFETY VALVE




In the event of a malfunction of the pressure switch and the pressure inside the tank exceeds 10,5 bar, the safety valve operates mechanically by opening for venting outside the compressed air and return the pressure to the

safety value.



UPON THE INTERVENTION OF THE SAFETY VALVE IT IS COMPULSORY TO TURN OFF THE MACHINE, NOTIFY THE SUPERVISOR AND REQUEST ASSI-STANCE OF MAINTENANCE TECHNICIANS OR AUTHORIZED DEALER TO MAKE THE CONTROL OF PRESSURE SWITCH.

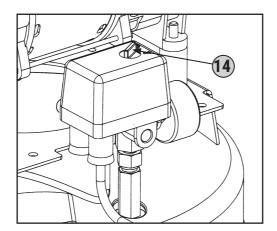



TROUBLESHOOTING



FAILURE TO FOLLOW THE PROCEDURES CONTAINED IN THIS PARAGRAPH MAY CAUSE A MALFUNCTION OF THE MACHINE AND/OR UNDERMINE THE SAFETY OF THE SAME.

PROBLEM	CAUSE	SOLUTIONS
The machine does not turn or stops and does not restart.	No electricity.	 Check that the switch (14) is in pos. "ON (I)". Check that the power cord plug (11) is working properly and inserted into the socket. Check that any extension cords used and the current outlet are operatio- nal. Check that the switch of the main po- wer supply is functioning and in pos. "ON (I)".
	Intervention of the re- settable thermal due to current overload and/or short circuit.	 Follow the procedure described in par. "RESETTABLE THERMAL".
The machine often starts again without the pneumatic plant is operating.	Air leaks from the pipe or its pneumatic plant.	 Check that the pneumatic system is not damaged. Check the integrity of the pipe con- nections and air system.

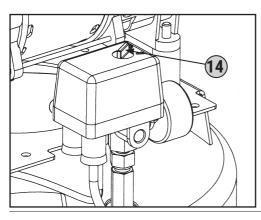

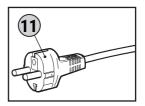

PROBLEM	CAUSE	SOLUTIONS
Intervention of the sa- fety valve.	Faulty pressure switch.	
Air leak from the valve of the pressure switch with machine stopped.	Check valve is dirty or worn.	See the maintenance technicians or contact an Authorized Reseller.
The machine vibrates and/or makes a lot of noise.		
Frequent starts and low efficiency.	Dirty air filters.	

MAINTENANCE



This is the set of activities designed to maintain the conditions of use and operation of the machine, through various types of intervention (adjustments, visual inspection, recovery levels, etc.) performed by the maintenance technician or by the manufacturer's technician at the frequency forecast.





THE AUTHORIZED OPERATORS SHOULD PERFORM ONLY THE MAINTENAN-**CE WORKS REQUIRED IN ACCORDANCE WITH THEIR SPECIFIC PROFESSIO-**NAL COMPETENCE AND UPON CONSENT OF THE SUPERVISOR IN CHARGE.



IT IS FORBIDDEN TO AUTHORIZED OPERATORS TO LEAVE THE MACHINE WI-THOUT SUPERVISION DURING THE MAINTENANCE WORKS. THEY SHOULD ALSO PROTECT THE PLACE OF WORK WITH A WHITE/RED CHAIN AND EXPO-SE ADEOUATE NOTICE OF DANGER AND PROHIBITION FOR WORKS IN PROGRESS.

	POINT OF INTERVENTION	TANK
	TYPE OF INTERVENTION	CONDENSATE DRAIN


EVERY WEEK



At the end of the workday, run the condensate drain of the tank, following the steps below:

- 1) Place under the condensate drain valve, a collection container fit for use.
- Slowly open the valve of condensate drain (9) and wait for a complete exit of compressed air and condensate.
- **3)** Close the drain valve.
- 4) Indicate the action on the maintenance card.

YOU MUST DISPOSE OF THE CONDEN-SATION LIQUID IN COMPLIANCE WITH THE REGULATIONS IN FORCE IN THE COUNTRY OF USE.

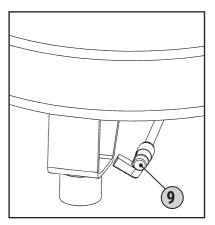

EVERY MONTH



Proceed as follows:

1) Unscrew the screw (1).

2) Remove the cover (2).

3) Remove the filter (3) and clean it with compressed air.




REPLACE THE FILTER IN CASE OF EVIDENT WEAR.

- 4) After placing back the filter, replace the cover properly.
- 5) Indicate the action on the maintenance card.



YOU MUST DISPOSE OF THE FILTERS REPLACED IN COMPLIANCE WITH THE REGULATIONS IN FORCE IN THE COUNTRY OF USE.

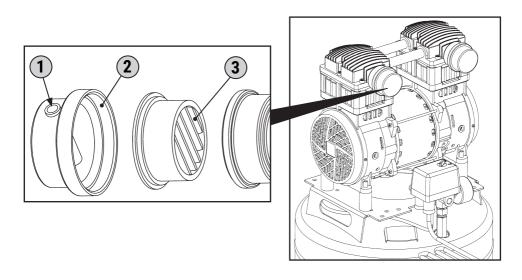


	POINT OF INTERVENTION	COALESCING FILTER
	TYPE OF INTERVENTION	REPLACE

DRYER FILTER ESC.F1

MAINTENANCE EVERY 1000 HOURS OR EVERY 12 MONTHS


- Replacement of PELLETS REGENERATING KIT (COD. 84213)
- Unscrew the dryer from the beam (17).
- Unscrew the nut (14) keeping the closing cap (2) locked.
- Place a box under the dryer to hold the pellets to be replaced.
- Remove the upper head (13) and extract the spiral (21), the outer tube (12) and the one internal (11).
- Empty the dryer completely from the exhausted pellets, collecting plates (9) and the spring (16).
- Unscrew the stay rod from the closing cap (2) and extract it from the lower head (5).
- Remove the closing plate (3), with the whole valve unit (24), from the lower head (5).
- Blow inside the pipes (11-12) and the heads (5-13) with compressed air.
- Replace all filters (6) and the O-Rings (1-4-7-8).
- Insert the closing cap (2) with the whole valve unit (24), in the lower head (5) with the closing plate (3) being careful to the position of the O-Rings (7-8).
- Reposition the tubes (11-12) in the lower head (5), insert the first plate (9) and fill the INTER-NAL tube (11) with new pellets, insert the second plate (9) and the spring (16).
- Insert the spiral (21) with the spacers (22) engaging it to the fitting (23).
- Reposition the upper head (13) being careful to the position of the threaded hole for the fixing it to the beam.
- Replace the copper washer (15) and tightening the nut (14) torque 25 Nm.
- Fasten everything to the spar (17) with screws (18), washers (19) and spacers (20).
- · Reconnect the system.
- Verify proper operation by doing attention to the pneumatic seal on the heads (5-13) on the external tube (12) and on the spiral couplings (21).

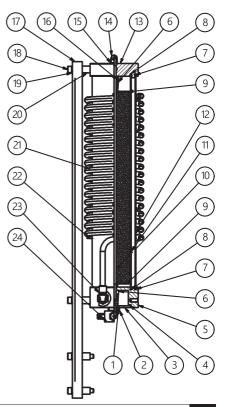

DRYER FILTER ESC.F3

MAINTENANCE EVERY 1000 HOURS OR EVERY 12 MONTHS

- Replacement of PELLETS REGENERATING KIT (COD. 84213)
- Unscrew the dryer from the beam (17).
- Unscrew the nut (14) keeping the closing cap (2) locked.
- Place a box under the dryer to hold the pellets to be replaced.
- Remove the upper head (13) and extract the outer tube (12) and then the one internal (11).
- Empty the dryer completely from the exhausted pellets, collecting plates (9) and the spring (16).
- Unscrew the stay rod from the closing cap (2) and extract it from the lower head (5).
- Remove the closing plate (3) from the lower head (5).
- Blow inside the pipes (11-12) and the heads (5-13) with compressed air.
- Replace all filters (6) and the O-Rings (1-4-7-8).
- Insert the closing cap (2) in the lower head (5) with the closing plate (3) being careful to the position of the O-Rings (7-8).
- Reposition the tubes (11-12) in the lower head (5), insert the first plate (9) and fill the INTER-NAL tube (11) with new pellets, insert the second plate (9) and the spring (16).
- Reposition the upper head (13) being careful to the position of the threaded hole for the fixing it to the beam.
- Replace the copper washer (15) and tightening the nut (14) torque 25 Nm.
- Fasten everything to the spar (17) with screws (18), washers (19) and spacers (20).
- · Reconnect the system.
- Verify proper operation by doing attention to the pneumatic seal on the heads (5-13) on the external tube (12).

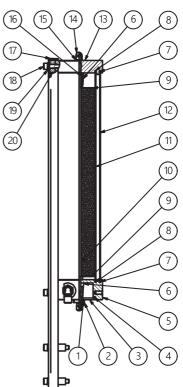



EXTRAORDINARY MAINTENANCE



This is the set of activities designed to maintain the conditions of use and operation of the machine, through various types of intervention (adjustments, changes, etc.) performed by qualified technicians in the event of failure or wear.





FOR ANY INTERVENTION OF EXTRAORDINARY MAINTENANCE IT IS COM-PULSORY TO REQUIRE THE TECHNICAL ASSISTANCE TO THE MANUFACTU-RER OR AUTHORIZED DEALER.

DECOMMISSIONING

STORAGE



In the event that the machine should be stored for some time, you must place it in a safe environment, with a proper degree of temperature, humidity and protected from dust.

DEMOLITION



🔊 Upon demolition of the machine, strictly follow the requirements in force in the country of demolition.

Proceed to the differentiation of the parts that make up the machine according to the different types of construction materials (plastics, copper, iron, etc. ...).



These products must be considered pollutants and hazardous charging licensed contractors and disposed of by appointing specialized companies for different types of products.



ALL PARTS OF THE MACHINE MUST BE IDENTIFIED ACCORDING TO THE DEFINITIONS OF "CER CODES" (EUROPEAN WASTE CATALOGUE) AND DI-SPOSED OF BY INSTRUCTING AUTHORIZED AND SPECIALIZED COMPANIES. IN STRICT COMPLIANCE WITH THE REGULATIONS IN FORCE IN THE COUNTRY OF

DEMOLITION OF THE MACHINE.



WASTES OF ELECTRICAL AND ELECTRONIC EQUIPMENT (RAEE) MUST BE **DISPOSED OF IN STRICT COMPLIANCE WITH THE REGULATIONS IN FORCE** IN THE COUNTRY FOR THE DEMOLITION OF THE MACHINE.

SPARE PARTS



GENUINE PARTS FOR ANY SUBSTITUTIONS MUST BE REQUIRED ONLY TO THE MANUFACTURER



IT IS STRICTLY PROHIBITED THE USE OF NON ORIGINAL SPARE PARTS.



Via delle Tezze, 20/22 36070 TRISSINO (VI) Italy

phone +39 0445 96 20 00 info@gentilinair.com

www.gentilinair.com



41011702 - 03/2022